

MFIS_2021 Installation notes

System requirements

OS Name Microsoft Windows 10
System Type x64-based PC
Processor Intel(R) Core(TM) i7 CPU 860 @ 2.80GHz, 2800 Mhz, 4 Core(s), 8 Logical Processor(s)
Installed Physical Memory (RAM) 32.00 GB (at least 16.00 GB is needed)

Download options


Please download MFIS software package:

[Seidel Software MFIS 2021.zip \(2.22 GB \)](#)

- AnI-3SF_12.07.2021
- Margarita (project) LV2020 x64SF
- Kristine_GLOBAL LV2020 x64_sf

Seidel_Software_ MFIS_2021

The downloaded *.zip file of the software should be directly extracted on to C:\ drive (important for proper functionality!!!): if you have software in the subfolder with the same name of *.zip file then move this folder directly on C:\ drive (see below).

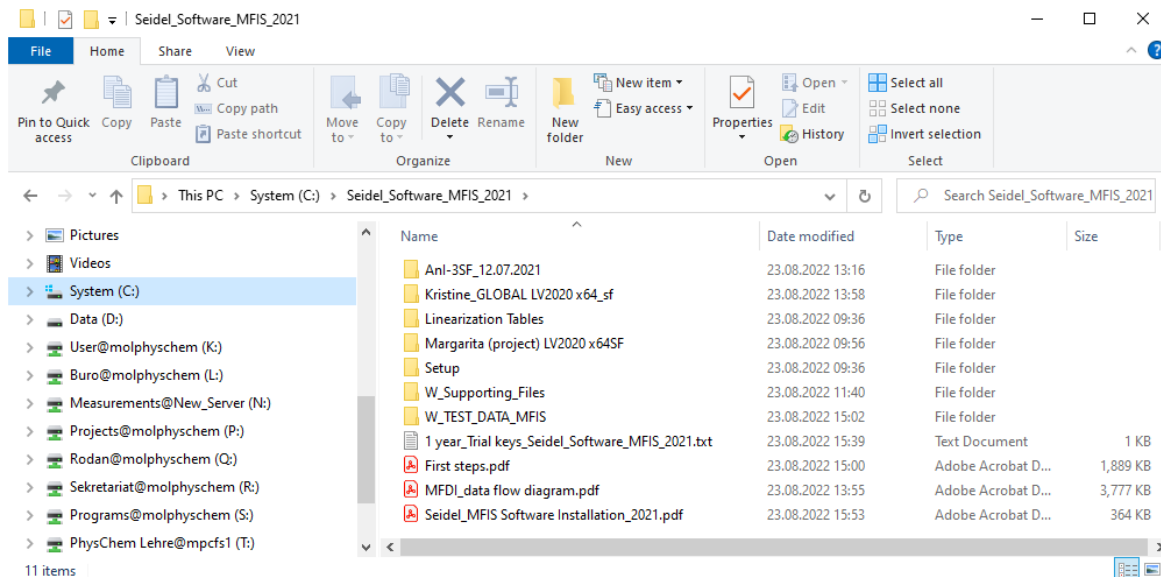

LV runtime Engines

- LabView2020x64 Runtime Engine (ni-labview-2020-runtime-engine_20.1.0_offline.iso) and LabView Vision RunTime (ni-vision-runtime_20.7.0_offline.iso, **National Instruments license is needed !!!**) should be installed by running corresponding setup files:

C:\Seidel_Software_MFIS_2021\W_Supporting_Files\LabView 2020 Runtime Engines

- The following files should be installed additionally for “AnI-3SF” software if you get error message about missing *.vi’s or *.dll’s by program start:

C:\Seidel_Software_MFIS_2021W_Supporting_Files\MASMTRunTimeEngine\ setup.exe

and

C:\Seidel_Software_MFIS_2021W_Supporting_Files\VC2017_redist.x64 \ VC2017_redist.x64.exe

and


C:\Seidel_Software_MFIS_2021W_Supporting_Files\NDP452-KB2901907-x86-x64-AIiOS-ENU\ NDP452-KB2901907-x86-x64-AIiOS-ENU.exe

and

C:\Seidel_Software_2017W_Supporting_Files\VC++ runtime 6195\ vcredist_x64.EXE

IMPORTANT!!!

US regional settings (decimal symbol - dot, digit grouping symbol - comma) are required to be set up on the computer.

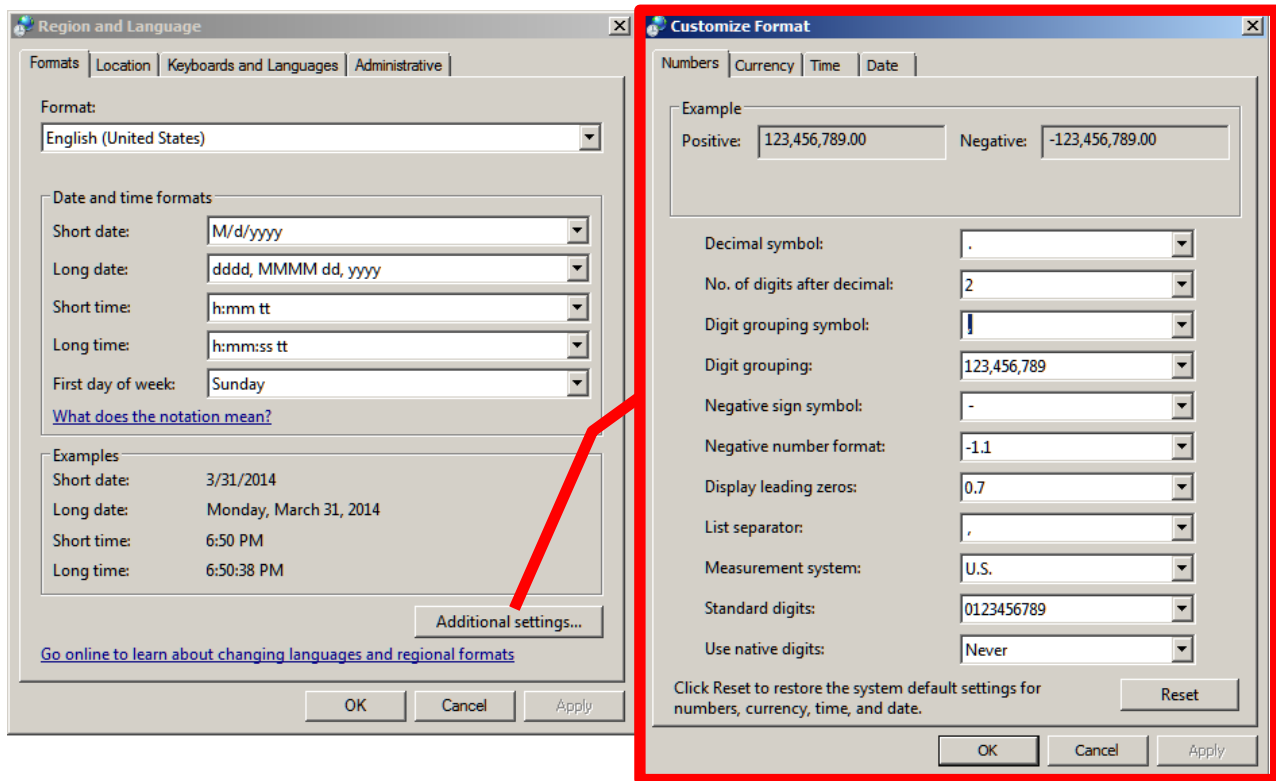


To run the software a trial license key (time limited) or permanent MAC key related to the physical address of PC (MAC) is needed!!!

For license key please send an email to Dr. Suren Felekyan: suren.felekyan@hhu.de

1 year valid trial keys are included in “1 year_Trial keys_Seidel_Software_MFIS_2021.txt” file !!!

For Margarita and Kristine software:

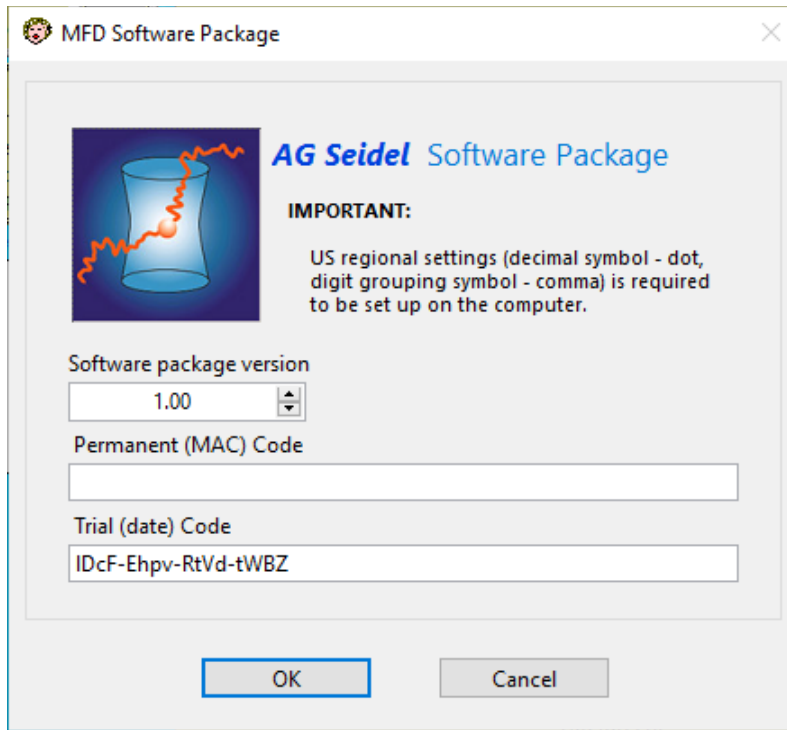

For Anl-3SF software:

